

Planmeca Compact™ i

Upgrade
Program



ENGLISH

PLANMECA

1. 2nd generation touch panel

Upgrade your **Planmeca Compact™ i** to the newest generation. This retrofit option includes a new powder-painted aluminum console frame with a fully integrated, state-of-the-art touch panel. Clear, large symbols and texts are displayed on the touch screen, making the use of your Planmeca Compact i even more efficient and enjoyable.

This kit also includes the multifunctional Planmeca Compact i foot control – more versatile than before.



2. Foot control

The foot control of **Planmeca Compact™ i** dental units provides users with great flexibility. All necessary functions can be carried out from it, keeping your hands free and minimising cross-contamination. The foot control is intuitive in use, and its significant weight ensures stable placement on the floor.

You can also use the foot control to operate your **Planmeca PlanScan®** intraoral scanner, if integrated with your dental unit.



3. Upholstery

Complement the splendid design of your dental unit by giving it a personal touch with your favourite colour. Select the perfect one from our inspiring collection and create the look of your dreams!



Our **Comfy™** upholstery is made of durable artificial leather. The thin structure ensures excellent reach to the treatment area, and the seamless design allows easy upkeep and enhanced infection control.

Made of extremely durable artificial leather and viscoelastic memory foam, our **Ultra Relax™** upholstery adapts perfectly to patients' body shapes, allowing them to stay relaxed even during the longest treatment sessions.

Comfy

Ultra Relax



Arm rests

4. Bowl

The dental unit bowl is subjected to frequent use in treatment situations. As it is heavily exposed to fluids and is one of the most visible parts of the unit, it is recommended to change the bowl assembly at times. With a different colour you can also provide your dental unit with a new look.



Bowl assembly, blue



Bowl assembly, clear sandblast



5. Balanced arm reach

Planmeca is known for excellent ergonomics. To ensure a better reach to instrument hoses, the pivot point of the balanced arms have been lowered. The reach of your instrument with the new balanced arm is approximately 4cm longer than with the previous arm. Furthermore, the arms are individually adjustable to compensate for the weight of the instrument. This is functionality at its best.



6. Headrest

The headrest mechanism is one of the most used parts of **Planmeca Compact™ i** dental units. The mechanism contains small fine-mechanical parts to enable effortless and correct positioning of the patient's head. We know from experience that the headrest will show wear and tear from frequent use. We are happy to offer you the possibility to replace the headrest mechanism for a brand new one with this top offer.



7. Planmeca Romexis® Clinic Management

Our **Planmeca Romexis® Clinic Management** software module is a solution no other manufacturer can offer. It allows Planmeca digital dental equipment to be easily connected to a network to generate valuable data and enable attractive new services for efficiency and quality assurance purposes. Smooth operation is guaranteed, as real-time information on dental units and their status is always readily available. This helps to locate and optimise available resources, and also provides the technical staff with the tools to thoroughly monitor the entire clinic with actual data.



8. Suction arm with Flexy™

The modular **Flexy™** holder for suction tubes and instruments supports all your treatment needs. You can also easily integrate your **Planmeca PlanScan®** intraoral scanner and tablet device for an optimal chairside workflow.

Adjustable suction arm
with Flexy holder



Chair mounted left/right suction arm
with Flexy holder



9. Infection control kit

Reliable infection control is of paramount importance to any dental unit. This extensive package will replace all critical parts and sealings to ensure maximum infection control throughout the treatment process. The kit includes filters, new suction tubes and suction tube handpieces.

10. Yearly maintenance kit

To guarantee optimal performance, it is crucial to conduct yearly maintenance on your dental unit. Our yearly maintenance kit includes e.g. necessary filters, valves, membranes, O-rings, and a new foot control cover.



A safer
treatment
environment

Always
clean
water


11. Suction Tube Cleaning System

Planmeca has designed an intelligent Suction Tube Cleaning System to guarantee proper functioning of the dental unit suction system, and ensure a safe treatment environment. The daily cleaning procedure is fast and simple due to the intuitive user-interface. The system comes with a new and enhanced selective suction system with double coarse filters. Cleaning your suction tubes has never been this simple or efficient.



12. Planmeca WEK


Planmeca WEK is an easy-to-use system that continuously feeds disinfectant into instrument and cup fill waterlines. The feeding is done automatically and users only need to replace the disinfectant bottle when it is getting empty. The system uses Metasys Green&Clean WK disinfectant.

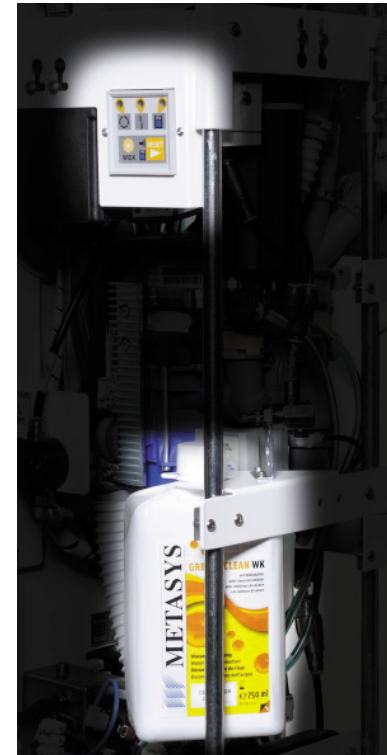

13. Planmeca SingLED®

An operating light convinces with its powerful illumination, and **Planmeca SingLED®** truly excels in this area. The operating light is powered by a single 5W LED that is reflected over a faceted reflector. This ensures an amazing shadow-free field of light with dazzle-free sharp edges. The three rotating axes provide a versatile movement range and the best possible positioning for each procedure. The large illumination field provides excellent visibility to every corner of the treatment area.

Planmeca SingLED is a must for any modern dental unit.

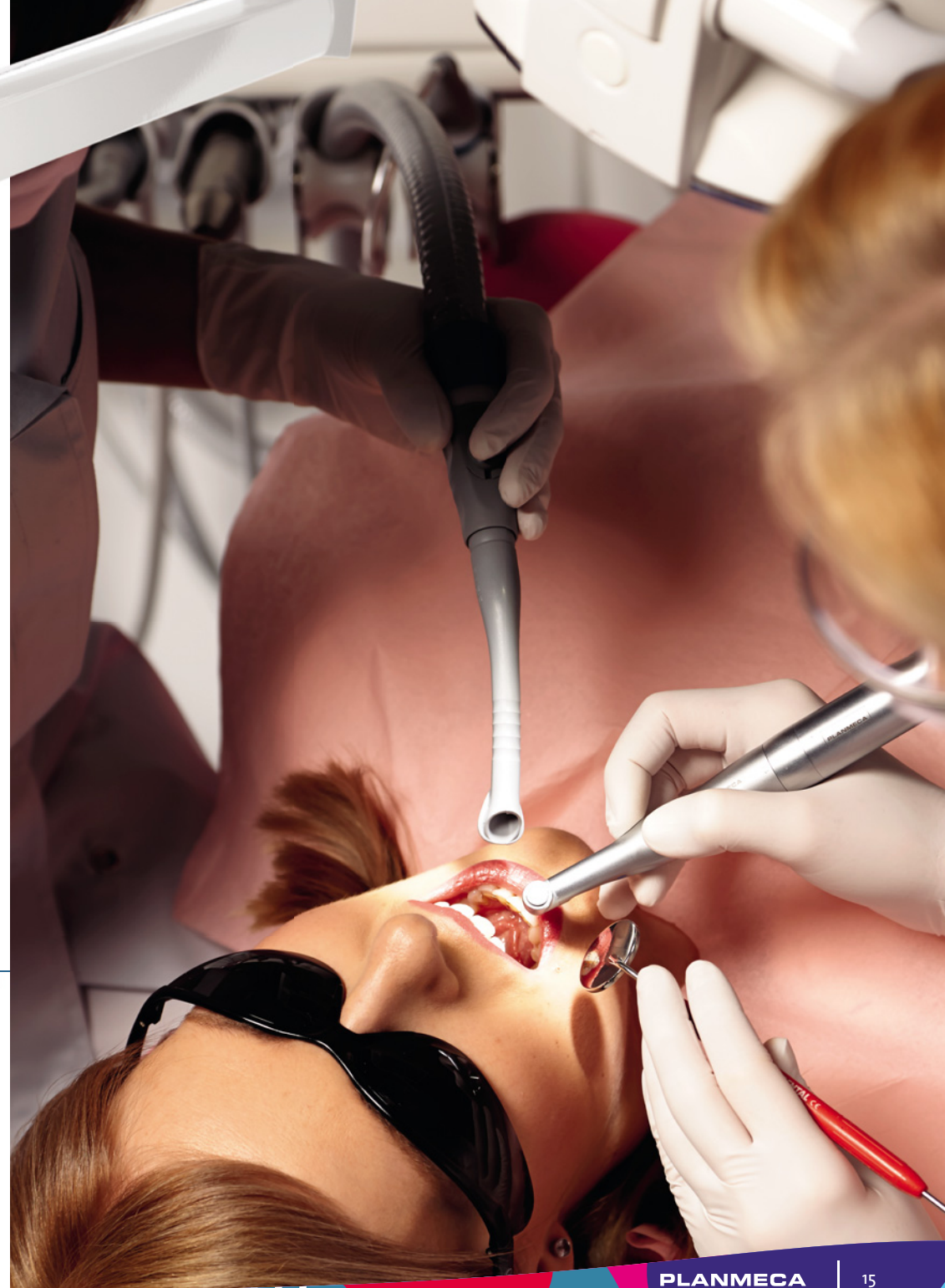



All retrofit options are only for Planmeca Compact™ i units manufactured from 2002 to 2011 (serial numbers starting with UPMC and UPIX). The dealer is responsible for informing Planmeca about any modifications done to the unit after manufacturing. All parts in these retrofit kits are non-refundable.

It is recommended that the dental unit software is always updated to the latest version when upgrading for new features. The latest software versions can be downloaded from Planmeca Dealer Support.

The campaign is valid until the end of 2016.

PLANMECA

Asentajankatu 6 | 00880 Helsinki | Finland | tel. +358 20 7795 500 | sales@planmeca.com | www.planmeca.com

Images may contain optional items not included in standard delivery. Available configurations and features may have country or area specific variations. Some products displayed above may not be available in all countries or areas. Rights for changes reserved.

Planmeca, All in one, Anomat Plus, Cobra, Comfy, DentoVac, Digital perfection, Economat Plus, Elegant, Flexy, Mini-dent, Perio Fresh, PlanEasyMill, Planmeca 4D, Planmeca AINO, Planmeca ARA, Planmeca CAD/CAM, Planmeca Chair, Planmeca Clarify, Planmeca Compact, Planmeca FIT, Planmeca Intra, Planmeca iRomexis, Planmeca Lumion, Planmeca Lumo, Planmeca Maximity, Planmeca Minea, Planmeca Minendo, Planmeca Minetto, Planmeca mRomexis, Planmeca Noma, Planmeca Online, Planmeca PlanCAD, Planmeca PlanCAM, Planmeca PlanClear, Planmeca PlanMill, Planmeca Planosil, Planmeca PlanPure, Planmeca PlanScan, Planmeca ProCeph, Planmeca ProFace, Planmeca ProID, Planmeca ProMax, Planmeca ProModel, Planmeca ProOne, Planmeca ProScanner, Planmeca ProSensor, Planmeca ProX, Planmeca Romexis, Planmeca Serenus, Planmeca SingLED, Planmeca Sovereign, Planmeca Ultra Low Dose, Planmeca Vision, Planmeca Verity, Planmeca Waterline Cleaning System, Planmeca Xtremity, Proline Dental Stool, ProTouch, Saddle Stool, SmartPan, SmartTouch, Trendy and Ultra Relax are registered or non-registered trademarks of Planmeca in various countries.

10039316/0316/en

